

Quill

BARBED
SUTURE

by corzamedical

The **exceedingly**
capable suture.

Eliminating suture knots is just the beginning of its capabilities.




Quill[®] knotless barbed suture eliminates the need to tie tedious and time-consuming knots.

Utilize a running stitch to make approximation faster with the strength of an interrupted stitch. The barbs lock the suture in place to maintain tension without assistance or recoil.

More anchoring surface area.

High Density (HD) barbs maximize tissue engagement and hold.

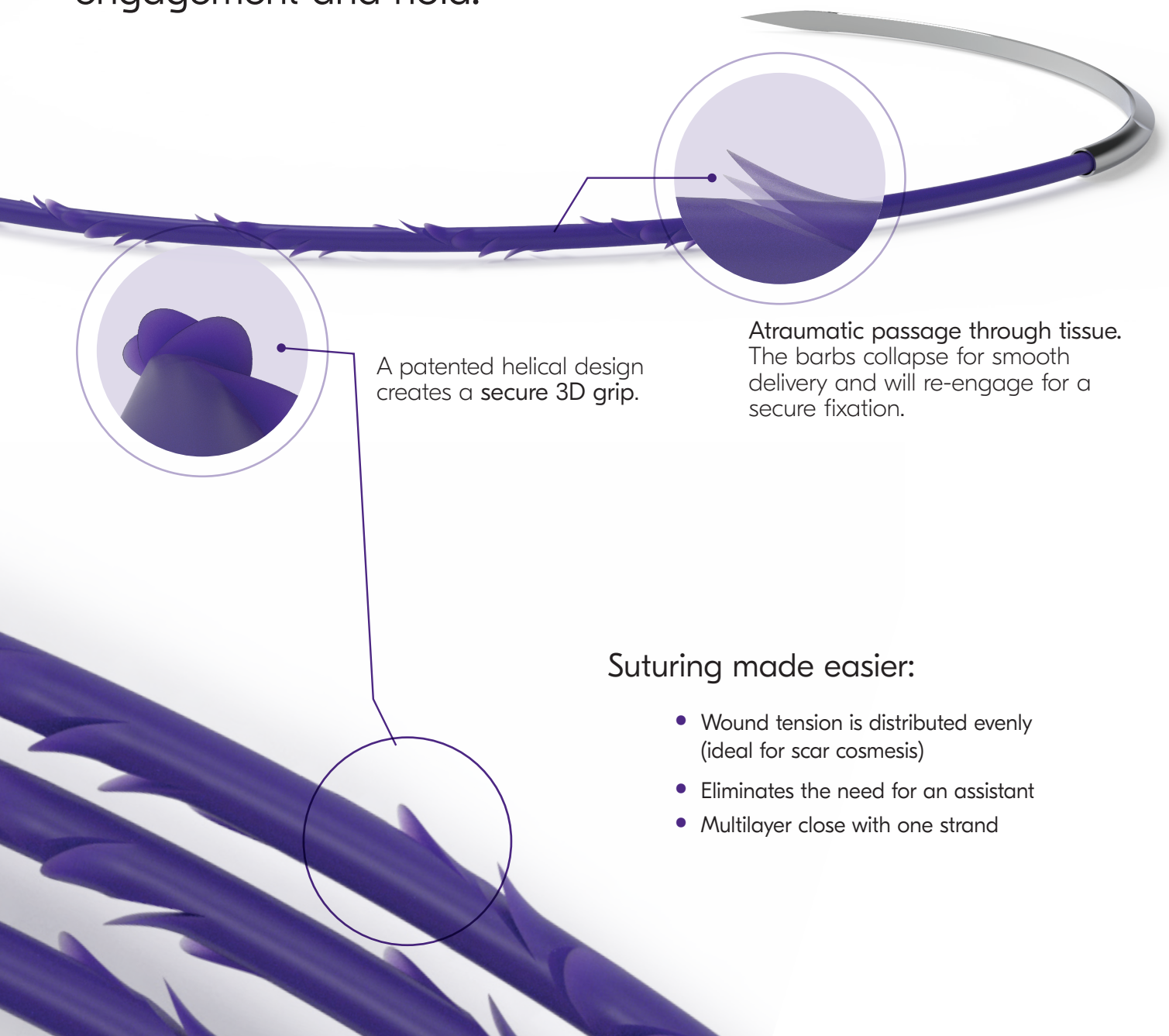

A patented helical design creates a secure 3D grip.

Atraumatic passage through tissue. The barbs collapse for smooth delivery and will re-engage for a secure fixation.


Suturing made easier:

- Wound tension is distributed evenly (ideal for scar cosmesis)
- Eliminates the need for an assistant
- Multilayer close with one strand

Quill barbed sutures offers a versatile portfolio.

Two configurations in a variety of lengths and USP sizes.

Bi-directional: Double-armed with a midpoint where barbs reverse direction.

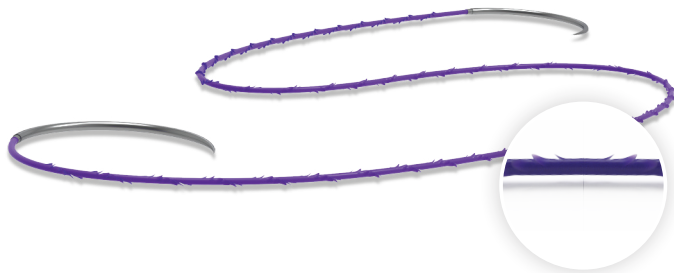

Bi-directional suture can be used simultaneously by two surgeons.

Uni-directional: Single armed with an adjustable anchoring loop.



The adjustable anchoring loop is large for the needle to easily pass through.



A wide array of needle options across multiple specialties, such as:

- Small reverse cutting needles for cosmetic and delicate surgeries.
- Small taper needles to maneuver in tight spaces for robotic surgeries.
- 5/8 circle needles specific for urologic procedures.

Suture Materials.

Short-Term Absorbable

Monoderm™ (PGA-PCL)

Undyed (Clear), Dyed (Violet)

Essentially complete in 90-120 days

Long-Term Absorbable

PDO (Polydioxanone)

Dyed (Violet)

Essentially complete in 180 days

Non-Absorbable

Polypropylene

Undyed (Clear), Dyed (Blue)

Permanent



Packaging made easier.




Color-coded foil blister for easy identification.

The single barrier foil blister reduces waste. Easy-to-read labels cover the entire pouch.

Minimal suture memory with the tray.

Environmentally friendly with recyclable polypropylene tray

Labeled tray for easy and quick identification in the sterile field

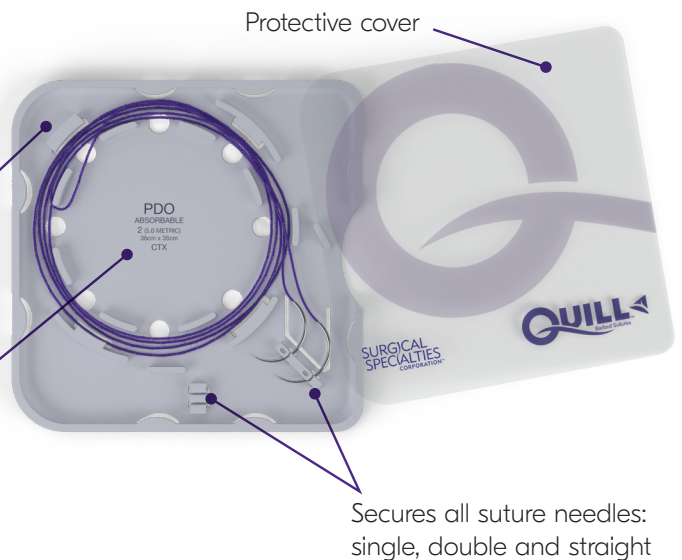

Start your free trial today.

service@corza.com

877-991-1110

www.corza.com

Quill is available through distribution or direct. For questions or a demonstration, please contact your local sales representative or contact us. For procedural videos, visit youtube.com/user/TheQuillDevice



Surgical Specialties, a Corza Medical company, 247 Station Drive, Suite NE1, Westwood, MA 02090, U.S.A.
All trademarks are property of their respective owners. © 2021 Corza Medical. All Rights Reserved. COR-220 R3 US 7/21